

STORAGE & PREPARATION GUIDE

AlloFuse[®] SELECT
CELLULAR BONE MATRIX

Natural **osteoconductive**, **osteoinductive** and **osteogenic** properties with **viable bone forming cells** to help initiate and nurture bone growth, delivering the benefits of an autograft without the drawbacks.

STORAGE

Upon receipt, remove AlloFuse[®] Select from its shipping container. Store the product in its original branded package at $-80^{\circ}\text{C} \pm 5^{\circ}\text{C}$ ($-112^{\circ}\text{F} \pm 9^{\circ}\text{F}$) to maintain product integrity and preserve its maximum 5-year shelf life.

AlloFuse[®] Select may be temporarily stored between -45°C to -75°C (-49°F to -103°F) for up to 90 days. To restore the full 5-year shelf life, return the product to storage at $-80^{\circ}\text{C} \pm 5^{\circ}\text{C}$ ($-112^{\circ}\text{F} \pm 9^{\circ}\text{F}$) at any time within the 90-day temporary storage period.



Ready to Use in a Few Simple Steps

PREPARATION

- Remove inner tray from the box. Using aseptic technique, peel open the tray from the marked corner.

- Transfer the vial into the sterile field. Place the closed vial in a sterile saline bath that is between 35°C (95°F) and 39°C (102°F).

NOTE: It is recommended to keep the vial closed during thawing. Open the vial just prior to use to ensure the bone allograft stays hydrated maintaining cell viability.

WARNING: Do not use a bath at a temperature greater than 39°C (102°F).

- Thaw the AlloFuse Select for approximately 10 minutes (1 cc) to 25 minutes (5 to 10 cc). Note that the 10 cc vial might take longer to completely thaw.


NOTE: For increased operating room efficiency, AlloFuse Select has up to 6 hours of working time after it is thawed. It is recommended to thaw AlloFuse Select during the beginning of the case.

NOTE: Use within 4 hours is required in New York.

- AlloFuse Select is now ready to use. Remove it from the vial and implant into place. Discard any unused product.

NOTE: Please refer to the AlloFuse Select package insert for more information.

AlloFuse® Select is ready for use after thawing – no decanting or mixing required.

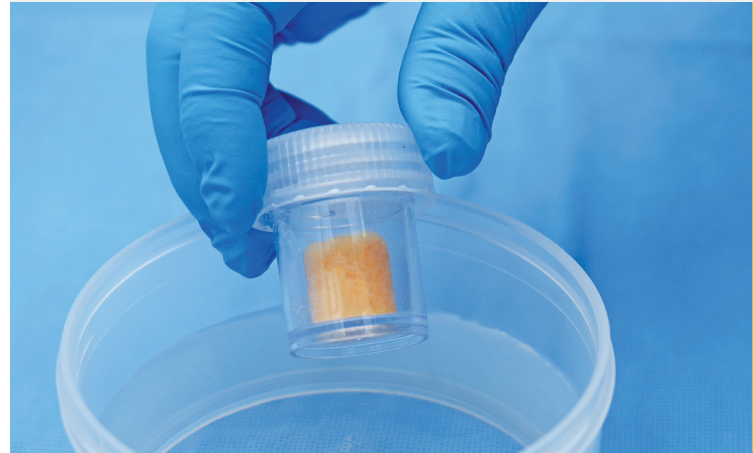

AlloFuse® SELECT

CELLULAR BONE MATRIX

ORDER NOW
800. 557. 3587

AlloSource, a life sciences organization, helps restore patient functionality by transforming the gift of human tissue donation into enhanced medical products that enable a life of movement, health, and wellbeing. We partner with medical professionals and biomedical companies who share our drive for innovation and quality in the use of human tissue allografts for medical advancements. Headquartered in Centennial, Colorado, we have served a global marketplace since 1994. Learn more at allosource.org.

6278 S Troy Cir
Centennial, CO 80111
USA

MAIN 720.873.0213
allosource.org

